


## SHIMADZU GC-2030/GC-2050

Clarity Control Module ENG

Code/Rev.: M240/100D Date: 2025-10-22

Phone: +420 251 013 400 clarity@dataapex.com www.dataapex.com DataApex Ltd.
Petrzilkova 2583/13
158 00 Prague 5
Czech Republic

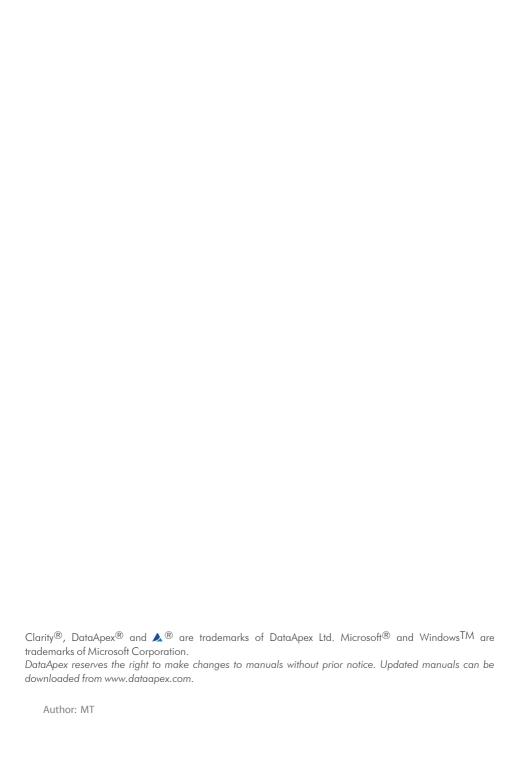

# **Contents**

1 Shimadzu GC-2030/GC-2050	1
2 Requirements	2
2.1 Software requirements	2
2.2 Hardware requirements	2
3 Installation procedure	
3.1 Installing Correct Version of ICF	3
3.2 Network connections	4
3.3 Installation Qualification of Shimadzu GC-2030/GC-2050	5
4 Clarity - Shimadzu Edition Configuration	8
5 Using the Shimadzu GC-2030/GC-2050	13
5.1 Method Setup - GC	
5.2 Device Monitor	17
6 Troubleshooting	18
6.1 Problems specific to the Shimadzu GC-2030/GC-2050	19

To facilitate the orientation in the **Shimadzu GC-2030/GC-2050** manual and **Clarity - Shimadzu Edition** chromatography station, different fonts are used throughout the manual. Meanings of these fonts are:

Open File (italics) describes the commands and names of fields in **Clarity - Shimadzu Edition**, parameters that can be entered into them or a window or dialog name.

WORK1 (capitals) indicates the name of the file and/or directory.

ACTIVE (capital italics) marks the state of the station or its part.

Chromatogram (blue underlined) marks clickable links referring to related chapters.

The bold text is sometimes also used for important parts of the text and the name of the **Clarity - Shimadzu Edition** station. Moreover, some sections are written in format other than normal text. These sections are formatted as follows:

Note:	Notifies the reader of relevant information.
	Warns the user of possibly dangerous or very important information.

### Marks the problem statement or trouble question.

Description: Presents more detailed information on the problem, describes its causes,

etc.

Solution: Marks the response to the question, presents a procedure how to remove it.

## 1 Shimadzu GC-2030/GC-2050

This manual describes the use of the **Shimadzu GC-2030/GC-2050** with the **Clarity - Shimadzu Edition** software.



Fig. 1: Shimadzu Nexis GC-2030

**Shimadzu** GC-2030/GC-2050 are controlled using RC.NET driver developed by **Shimadzu** running under the Agilent ICF (Instrument Control Framework) library developed by Agilent, which can be run within the **Clarity - Shimadzu Edition** environment.

# 2 Requirements

## 2.1 Software requirements

**Shimadzu GC-2030/GC-2050** RC.net driver is exclusively available in **Clarity - Shimadzu Edition OEM** version. This OEM version is only available in selected countries in Europe region and sold through approved distributors. In terms of **Clarity - Shimadzu Edition** product numbers, this includes:

- Clarity Shimadzu Edition (p/n C50-28)
- GC Control license (p/n A23)
- GC Shimadzu license (p/n A23-001)
- optionally AS Control License (p/n A26) when the system contains an autosampler

Caution

A separate GC Shimadzu Control license (p/n A23-001) must be purchased for each controlled Shimadzu GC-2030/GC-2050 system. A similar license is required for **Shimadzu Nexera LC** systems.

**Shimadzu GC-2030/GC-2050** requires **Microsoft .NET version 4.8** or higher for correct installation and operation. The version is already installed on majority of PCs.

For complete list of .NET requirements, see the .NET Framework System Requirements on Microsoft web page.

Supported operating systems:

- Windows 10 (32/64 bit)
- Windows 11

Note:

Before installing Clarity - Shimadzu Edition, ensure that your Windows OS is updated to the latest version for optimal compatibility.

## 2.2 Hardware requirements

- A PC with a LAN interface is required.
- To ensure maximum compatibility, the module should run the latest firmware version (minimum 1.09).

For complete list of supported Windows OS and instrument firmwares, see Shimadzu webpages.

# 3 Installation procedure

**Shimadzu GC-2030/GC-2050** is part of the **Clarity - Shimadzu Edition** *Typical* installation. To install it, in the *Choose Components*select the *Typical* installation or check the Shimadzu GC-2030/GC-2050 in *Instrument Control Framework (ICF)* section during the installation of **Clarity - Shimadzu Edition**.

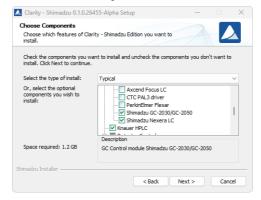


Fig. 2: Typical installation of Clarity - Shimadzu Edition

## 3.1 Installing Correct Version of ICF

Clarity - Shimadzu Edition expects a specific version of Agilent ICF. Because other programs may also be using Agilent ICF, it is possible they've installed a different version than is supported by Clarity - Shimadzu Edition. In that case in the installation you will be prompted to reinstall it in order to install correct version. Going forward with this step is crucial for correct functionality of Clarity - Shimadzu Edition and Shimadzu GC-2030/GC-2050.

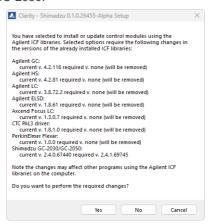


Fig. 3: Installing correct version of Agilent ICF

Note:

This reinstallation of Agilent ICF may cause that other programs using it, may not function properly.

### 3.2 Network connections

The instrument Shimadzu GC-2030/GC-2050 has to be connected to a site network by LAN. It is recommended to attach the hardware directly to the PC avoiding hubs, switches etc. Always contact your local LAN administrator who can make the appropriate settings.

Caution:

**Cross LAN** cable is primarily used for the direct connection of the instrument and the PC. This cable can also be used for the connection of the device to the switch or network socket, but with older switches, the **straight LAN** cable might be necessary.

### **LAN Settings**

PC: LAN card, TCP/IP protocol.

Both PC and Shimadzu GC-2030/GC-2050 should be configured on the same IP range.

### Firewall

Ensure that the firewall does not block communication from the Shimadzu GC-2030/GC-2050.

### 3.3 Installation Qualification of Shimadzu GC-2030/GC-2050

Shimadzu GC-2030/GC-2050is controlled via **Agilent ICF**, an external program developed by Agilent, and for that reason it must be validated using their utility. If you have installed Clarity - Shimadzu Edition with Agilent ICF, **the IQ is valid only if successful validation of ICF is attached**.

The validation of ICF can be performed directly from the IQ Report.

Caution:

When Clarity - Shimadzu Edition expects *Agilent ICF* is installed then **IQ** expects the same. If (due to any reason) the *Agilent ICF* installation is not found within Clarity - Shimadzu Edition, the *Installation Qualification Test* status is set by default to *FAILED*. To resolve this, it is necessary to re-install **Agilent ICF** through Clarity - Shimadzu Edition reinstallation and then perform **IQ** again.



Fig. 4: IQ Report with ICF installation present

Click the link "here" ①, after that it is necessary to click Run in two pop-up windows. Agilent Software Verification Tool window will open. Select what report file type should be generated and define post-qualification actions. The Reports folder field shows the location where the reports will be saved (C:\SVReports\). The folder can be opened directly by clicking the Open in the Agilent Software Verification Tool window. Click Qualify ② to run the IQ.

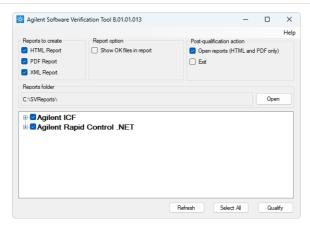


Fig. 5: Agilent Software Verification Tool

The HTML and PDF reports are opened automatically in the default browser/program if the *Open reports* option was enabled.

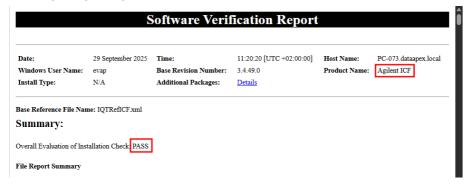


Fig. 6: ICF Report - PASS

At the end of the generated report, a table of installed drivers and their versions is displayed.

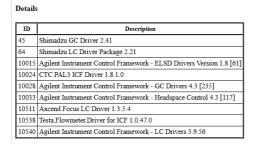


Fig. 7: ICF Report - Installed drivers table

The individual components can also be found in the component list if the *Show OK files in report* checkbox is selected in the *Agilent Software Verification Tool* window.

# 4 Clarity - Shimadzu Edition Configuration

- Start the Clarity Shimadzu Edition station by clicking on the licon on the desktop.
- Open the System Configuration dialog accessible from the Clarity Shimadzu Edition window using the System - Configuration... command.

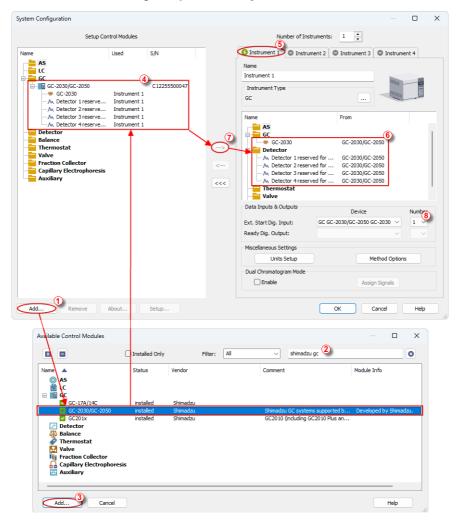


Fig. 8: System Configuration

- Press the Add button (See 4 on pg. 8.) to invoke the Available Control Modules dialog.
- You can specify the searching filter ② to simplify the finding of the driver.
- Select the **Shimadzu GC-2030/GC-2050** item and press the *Add* ③ button. The *GC-2030/GC-2050 Setup* dialog will open.

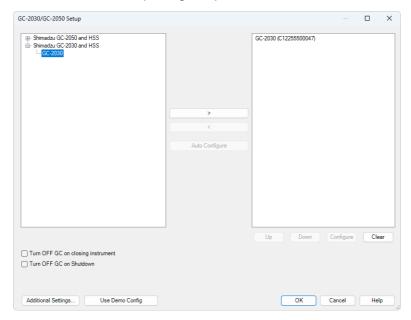


Fig. 9: GC-2030/GC-2050 Setup

Note: The Turn OFF GC on closing instrument and Turn OFF GC on Shutdown checkboxes govern the behavior of the GC while the Instrument window is closed and when Shutdown event is invoked (by user, by error or as a reaction from Event Table). When the function is not checked Clarity - Shimadzu Edition will not

be able to shut down the Shimadzu GC-2030/GC-2050.

• Select appropriate instrument type (a) and use Auto Configure (b) button to open *Automatic Configuration parameters* dialog.

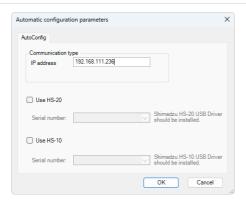


Fig. 10: Automatic configuration parameters dialog

- Insert the *IP address* of Shimadzu GC-2030/GC-2050 and confirm. The address can be found or set using instrument's touch panel.
- Then click *OK* button in the lower part of the *GC-2030/GC-2050 Setup* dialog and the instrument will then appear in *System Configuration*, including all modules of the instrument (detector, autosampler etc.).
- After previous steps drag and drop the device icon from the Setup Control Modules list ① on the left side of the System Configuration dialog to the desired Instrument ⑤ tab on the right side ⑥ (or use the -> button ⑦ to do so).
- Ext. Start Dig. Input setting <sup>®</sup> determines how Clarity handles external start signals from connected instruments:
  - "--" (dash) option Use this if your autosampler and detectors are interconnected by synchronization cables (or built-in together) so that the information about performed injection is passed on without direct input from Clarity Shimadzu Edition. Clarity Shimadzu Edition will initiate recording when it detects that the detector has begun acquiring data. This is the most common setup with Shimadzu GC-2030/GC-2050.
  - "1" option Use this setting if your autosampler and detectors are not directly interconnected by cables. Clarity Shimadzu Edition will coordinate the measurement by receiving a signal when the autosampler finishes injection and then triggering the detectors to start measuring. This applies mainly if your GC has an autosampler detectors are not built-in, so it is not so common.
- More detailed configuration parameters can be opened by using Configure button in GC-2030/GC-2050 Setup. In Configure GC-2030 window double click on specific module. More settings can be set in opened window, specific to the module.



Fig. 11: Configure GC-2030 dialog

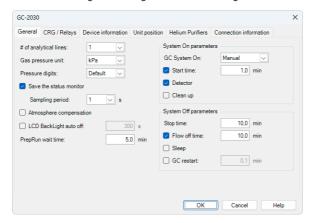


Fig. 12: Module configuration dialog

Note: Press the F1 key to display the original Shimadzu GC help with detailed


description of the dialog.

Note: The configuration dialog of the Shimadzu GC-2030/GC-2050 subcomponents

(detector/s, inlet...) can be displayed any time by double-clicking on its icon.

### Additional Settings...

This dialog allows you to configure additional settings of the control module.

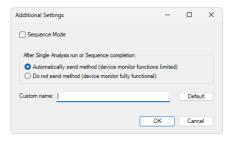

### Sequence Mode

Allows the autosampler to prepare the next vial while the previous run is still ongoing. Clarity - Shimadzu Edition will automatically lock one row ahead in the sequence.

To use this feature, the *Enable Sample Overlap* setting must be enabled in the ICF GC method (GC tab  $\rightarrow$  ALS section  $\rightarrow$  Tray/Other).

### After Single Analysis run or Sequence completion

Choose how Clarity - Shimadzu Edition sends the method to the device:

- Automatically send method Clarity Shimadzu Edition sends the method
  after each single analysis run or after the sequence is finished automatically.
  Device Monitor functions after the method is sent may be limited in this
  mode. This mode is necessary if you want to run a continuous loop of Single
  Analysis measurements. It does not affect Sequence runs.
- Do not send method The method is not sent automatically, allowing full access to Device Monitor functions even when method is sent to the device.

#### Custom name

Allows you to rename the module. The custom name will appear in the *System Configuration* dialog including the *Data Inputs & Outputs* section.

If left blank, the module will use its default name. Use the Default button to reset the name at any time.

## 5 Using the Shimadzu GC-2030/GC-2050

This chapter goes through the setting of the parameters of the **Shimadzu GC-2030/GC-2050** in the **Clarity - Shimadzu Edition**.

Caution:

Before opening the *Instrument* window with configured **Shimadzu GC-2030/GC-2050** device, ensure there is not any other PC connected to GC. Otherwise there will be communication error upon opening.

## 5.1 Method Setup - GC

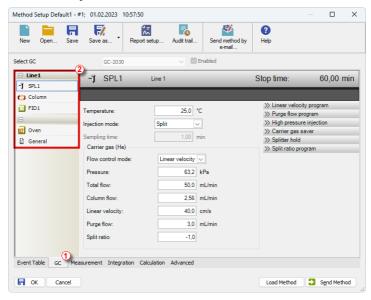


Fig. 13: Method Setup - GC

Note:

Press the *F1* key to display the **Shimadzu GC-2030/GC-2050 help** with detailed description of the dialog.

From **GC tab** ① you can set all parameters for individual subcomponents (detector/s, oven, inlet...) of instrument configuration.

For advanced setting of all other parameters use all tabs ② on the left side of the window in GC tab.

If there is configured any autosampler of AOC-20 family its parameters are about to set on **GC tab** in its respective section.

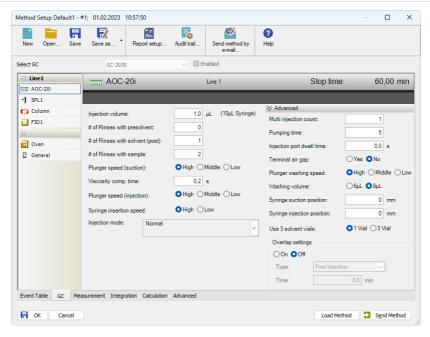


Fig. 14: Method Setup - GC - AOC autosampler

Note: Actual parameters read out from GC are used for new method or method adaptation instead of default parameters.

Note: For common analysis it is recommended to relate Oven temperature program via item Link oven program and acquisition time otherwise duration of the acquisition is governed by the item Stop time in the detector section.

If Carrier Gas Saver item is set On in the used method it is consequently required to trigger Start Prep Run manually by pushing Prep Run button above display of the GC or by some external contact. Any of these will transfer the GC to Ready state which is has to be reached prior actual analysis. If some external contact is about to be used the Start Prep Run has to be set as External in General section of the method. Besides that, it is necessary to have correctly set Relays in the "Module configuration dialog" on page 11 in the System Configuration. Refer to images below.

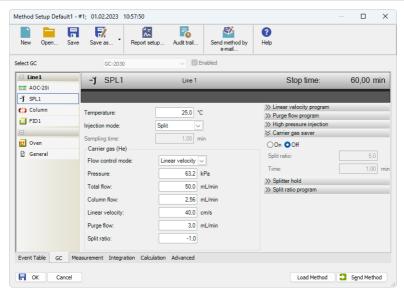


Fig. 15: Method Setup - GC - Gas Saver

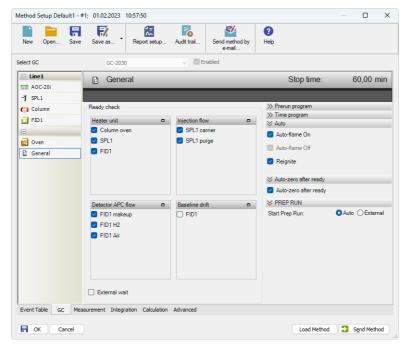


Fig. 16: Method Setup - GC - External Prep Run

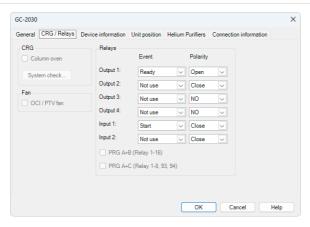


Fig. 17: Module configuration - CRG Relays

### 5.2 Device Monitor

The Device Monitor window can be opened by the Device Monitor command from the Analysis menu or using the **Device Monitor** @ icon in the *Instrument* window. You may use (icon for accessing Device Monitor from all Clarity - Shimadzu Edition windows. In the Device Monitor window, it is possible to monitor the status of the connected instruments.

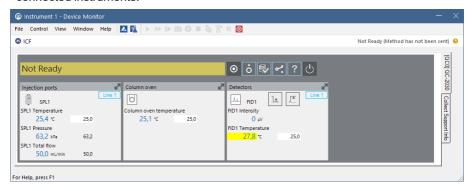


Fig. 18: Device Monitor

Press the F1 key to display the Shimadzu GC-2030/GC-2050 help with detailed Note: description of the dialog.

Note: Tabs on right side of the Device Monitor allowing setting collection of diagnostic information are displayed only when Clarity - Shimadzu Edition is started under

Administrator Windows account

Note: The Shutdown command works as expected only if item Turn OFF GC on Shutdown in the "GC-2030/GC-2050 Setup" on page 9 dialog of the System Configuration is checked. If it is checked, and Shutdown is performed, the setting defined in section System Off parameters from the tab General of the "Module configuration dialog" on page 11 dialog will be applied.

Some errors cannot be cleared using Device Monitor and have to be cleared via Note: manual commands using display of the GC.

- 17 -

# 6 Troubleshooting

When the solution for some problem cannot be discovered easily, the recording of communication between **Clarity - Shimadzu Edition** and Shimadzu GC-2030/GC-2050 control module can significantly help the **DataApex** support to discover the cause of the problem. The created \*.TXT files will greatly help in diagnosis of unrecognized errors and problems.

The recording can be enabled by adding or amending the LOGGING.INI file in the Clarity - Shimadzu Edition installation directory (C:\CLARITYSHIMADZU\CFG by default). The file can be edited in any text editor (e.g. Notepad). Following section should be edited:

[Log]
echo = ON
filename = log\_%D.txt
reset = OFF
; Sections List:
AuditTrail = ON
BadTrace = ON
CommandLine = ON
Acquisit = ON
AgilentICF = ON
SST = OFF
Internet = OFF
FractionCollector = OFF

Note:

%D (or %d) in the filename parameter means that the log will be created separately for each day. The *reset = OFF* parameter disables deleting the content of the log each time the station is started during the same day.

In case you cannot establish communication with instruments connected through ICF, please review the following issues:

### Check the network connection using the Ping command

The problem in communication between **Clarity - Shimadzu Edition** and Shimadzu GC-2030/GC-2050 may be caused by wrong network configuration, firewall preventing the connection, etc. Run the command line in Windows (for example by pressing the **Windows key** together with the **R** key, in the displayed *Run* window type *cmd* and press *Enter*).

In the command line type ping <ip-address-of-instrument> and press *Enter*. The *IP Address* is the same you entered during the configuration.

## 6.1 Problems specific to the Shimadzu GC-2030/GC-2050

Clarity - Shimadzu Edition can't be run and it displays "Agilent ICF is not installed correctly." message.

Cause: The cause of the problem is that the Agilent ICF has a different version than

expected by Clarity - Shimadzu Edition. It can typically happen when other software also using Agilent ICF decides to reinstall it. Thus next time Clarity -

Shimadzu Edition expects different version than is installed.

Solution: Solution is to reinstall Agilent ICF during Clarity - Shimadzu Edition installation.

Follow steps described in the chapter "Installation procedure" on pg. 3.

Print of Injection Control is not functional.

Cause: Agilent ICF is incompatible with some Clarity - Shimadzu Edition printing

procedures.

Solution: Method parameters of autosampler control are printed together with another GC

control method parameters using Instrument Control item within Method section in

Report Setup dialog.

ICF problems during installation or operation.

Cause: The cause of the problem might be that Microsoft .NET Framework is not enabled.

Agilent ICF requires Microsoft .NET Framework enabled for its function.

Solution: Check if Microsoft .NET Framework is enabled in Turn Windows features on or off dialog. Turn Windows features on or off dialog is accessible in Control Panel window under section Programs in its subsection Programs and Features. The

Microsoft .NET Framework related items are disabled by default in some Windows versions.

Programs and Features X III « Programs » Programs and Features Search Programs and Features ۵ Control Panel Home Uninstall or change a program View installed updates To uninstall a program, select it from the list and then click Uninstall, Change, or Repair. Turn Windows features on or Organize 🕶 ■ -Install a program from the Name Publisher Installed On Size network

Fig. 19: Location of Turn Windows features on or off in Windows 10

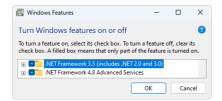


Fig. 20: Turn Windows features on or off - Windows 10